

Slings and lifting accessories

Guldmann™

Slings and other lifting accessories are a key part of the extensive range of Guldmann equipment for lifting, moving and positioning people.

Specially designed slings make it easy, practical and quick for the carer to support the user while being lifted and moved from one point to another, as well as ensuring that the carer can undertake such operations with a minimum of effort and strain.

All Guldmann slings are designed for use with both mobile and ceiling-mounted hoist systems, and are ideal for moving people to and from beds and wheelchairs, as well as in sitting and lying positions, as required.

Special Guldmann sling designs are available for lifting and moving users with

- different abilities to participate actively
- different degrees of body control and muscular strength
- different degrees of mobility
- single and double amputations
- special needs

and for both adults and children of many different sizes and weights. They can also be used as integral parts of training and activity programmes for users. All Guldmann slings are easily washable, to help ensure the best possible hygiene.

Accessories

In addition to slings, Guldmann provides a wide range of specially designed lifting accessories for helping both the user and carer ensure greater mobility with less effort. At the same time, they allow the carer to provide the best possible assistance, while working in ergonomically correct positions.

taylordolman

tel: 0845 838 0600

email: info@taylordolman.com



ABC slings

Guldmann™



1. Coloured loops

The coloured loops ensure symmetry in the lift

2. Material

Dirt repellent polyester characterized by its strength and elasticity

3. Padding on legs

All slings with divided leg supports have inset foam pads on legs. Easier to put on – smoother on skin

4. Label for registration of inspections

5. Label with size indication

Makes it easy and simple to find the correct size of sling

6. Instruction label

Indicates all that the carer needs to know

7. Padding in back section

For added comfort on models with small back surface

8. Auxiliary straps

It is possible to use auxiliary straps

9. Type marking

The colour of the band identifies type of sling

10. Positioning pocket

All Basic slings are provided with a pocket to ease correct positioning of sling

Fit

Improved fit and more sizes to choose from for added comfort and safety

Strength

Max. capacity increased to 250/450 kg

Washing

Easy handling, the sling can be washed at 85°C and tumble-dried

ABC slings Active

Guldmann™

1



2



1

Active Trainer

Function

Provides support around the chest and shoulder areas. Detachable leg straps give support around the hip area and prevent the user from sliding out of the sling.

Suitable for

- persons with reduced body balance, but must be able to bear weight on their legs
- walking training • supported standing and upright postures.

..... 2830xx

2

Active Micro Plus

Function

Supports the upper part of the body (from upper pelvis till under the shoulders) around the chest and upper thigh. The belt has a re-tightening function if necessary to adjust the support around the body.

Suitable for

- lifting persons with head control, but with reduced muscular strength of upper part of the body and above hip joint and thighbone
- in situations where the sling should be put on and taken off easily
- lifting to and from seated position
- in certain situations also lifting to and from a half-lying position together with Guldmann head support
- situations where a large opening at the bottom is required, such as access for toileting and hygiene.

Kids 2810xx

XS – XXL 2815xx

ABC slings Basic

Guldmann™

3



4



5



3

Basic Basic

Function

Supports the whole body excluding the head.

Suitable for

• lifting persons with head control and with or without body control • lifting to and from a seated or half-seated position (bed and wheelchair etc) • lifting and transfers from all patient furniture surfaces.

..... 2700xx

4

Basic Low

Function

Supports the body up to shoulder height.

Suitable for

• lifting persons with head and body control • lifting to and from a seated position, (chair and bed etc.) • lifting and transfers from all patient furniture surfaces.

..... 2710xx

5

Basic High

Function

Supports the whole body including the head.

Suitable for

• lifting persons with little control over their head and body • lifting to and from a lying position (bed and floor etc.) • lifting to and from a seated position • lifting and transfers from all patient furniture surfaces and the floor.

..... 2720xx

ABC slings Custom

Guldmann™

6



7



8



6

Custom Sit-On

Function

Supports the thighbone and back. The user sits on the sling. Added bands fixed across the opening at each side of the sling prevent the user from falling out sideways. Not suitable for persons with pressure sore problems.

Suitable for

- persons who remain on the sling for long periods due to high weight or customized wheelchair.

Other features

- helps reduce the risk of sliding out of wheelchair • material prevents heat generation between body and sling • sling is fitted on the person while lying, but is used when lifting to and from a seated position in bed or wheelchair.

..... 2910xx

7

Custom Sit-On High

Function

Supports the thighbone, back and head. The user sits on the sling. Added bands fixed across the opening at each side of the sling prevent the user from falling out sideways. Not suitable for persons with pressure sore problems.

Suitable for

- persons who remain on the sling for long periods due to high weight or customized wheelchair.

Other features

- helps reduce the risk of sliding out of wheelchair • material prevents heat generation between body and sling • sling is fitted on the person while lying, and is used when lifting to and from a lying position.

..... 2920xx

8

Custom Amputee

Function

Supports the whole body up to shoulder height and the pelvis region.

Suitable for

- lifting both single and double (above knee joint) amputees • lifting double and high amputees • lifting from a lying position. Head support should be used to support head.

..... 2900xx

ABC slings Kids

Guldmann™

9



10



9

Basic Basic

Function

Supports the whole body excluding the head. The inside belt offers an extra support of the body.

Suitable for

- lifting persons with head control and with or without body control
- lifting to and from a seated or half-seated position (bed and wheelchair etc.).

..... 2700xx

10

Basic High

Function

Supports the whole body including the head. The inside belt offers an extra support of the body.

Suitable for

- lifting persons with little control over their head and body
- lifting to and from a lying position (bed, changing bench and floor etc.)
- lifting to and from a seated position.

..... 2720xx

ABC Slings

Guldmann™

Size		Kids	Kids	Kids	XS	S	M	L	XL	XXL	XXXL
		4 - 6	6 - 10	10 - 14							
ACTIVE											
Micro plus	polyester	281001	281011	281021							
Micro plus	polyester				281531	281541	281551	281561	281571	281581	
Trainer	polyester	283001	283011	283021	283031	283041	283051	283061	283071	283081	
BASIC											
Basic	polyester	270001	270011	270021	270031	270041	270051	270061	270071	270081	
Basic	net	270102	270112	270122	270132	270142	270152	270162	270172	270182	
Basic Low	polyester				271031	271041	271051	271061	271071	271081	
Basic Low	net				271132	271142	271152	271162	271172	271182	
Basic High	polyester	272001	272011	272021	272031	272041	272051	272061	272071	272081	
Basic High ^{x)}	polyester										273031
Basic High	net	272102	272112	272122	272132	272142	272152	272162	272172	272182	
CUSTOM											
Amputee	polyester					290041	290051	290061			
Sit-On ^{xx)}	net						291052			291082	
Sit-On High ^{xx)}	net						292052			292082	

Max lifting capacity for all Guldmann ABC slings: 250 kg


x) Basic High XXXL can lift up to 450 kg

xx) Sit-On slings are available in two all-size models: Standard and Wide

CE

Guidelines on choosing the size of a sling

	Kids	Kids	Kids	XS	S	M	L	XL	XXL	XXXL
	4 - 6	6 - 10	10 - 14							
H	45 - 60	55 - 70	68 - 80	75 - 82	80 - 87	85 - 92	90 - 97	95 - 102	95 - 102	95 - 102
W	25 - 30	28 - 35	33 - 38	33 - 37	36 - 40	39 - 43	42 - 45	45 - 49	48 - 52	54 - 58

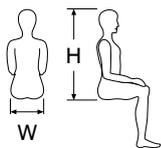

Materials




Nylon



Special woven polyester



Polyester

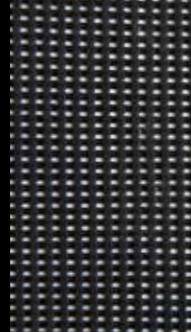

Net



The coloured loops ensure symmetry in the lift and at the same time ease the communication in connection with the lift.

Uni-D slings

Guldmann™

1



2



3



1



2



1

Uni-D

Function

Supports the whole body, excluding the head. This is a single patient use sling that tolerates water. Dispose after discharge. Beneficial for infection control program. Same design as Guldmann ABC Basic Basic sling.

Suitable for

- lifting persons with head control, but with reduced body control
 - lifting to and from a seated or half-seated position (chair and bed etc.).
- Packed in boxes of 10 slings.*

Material

Nylon

M/L	28362
L/XL	28363

2

Uni-D High

Function

Supports the whole body including the head. This is a single patient use sling that tolerates water. Dispose after discharge. Beneficial for infection control program. Same design as Guldmann ABC Basic High sling.

Suitable for

- lifting persons without head and body control
 - lifting to and from a lying position (bed, floor etc.)
 - lifting to and from a seated position.
- Packed in boxes of 10 slings.*

Material

Nylon

M/L	28364
L/XL	28365

3

Uni-D Turner

Function

To be used when turning a user into lateral position. This is a single patient use sling that tolerates water. Enables placing of sling under the patient without having to roll or move them. Smaller surface area of material for this sling allows access to patient's skin for inspection, treatment and care. Same design as Guldmann Turner. *Packed in boxes of 10 slings.*

Material

Nylon

..... 287005

Disposable slings

Guldmann™



1

Disposable High

Function

Supports the whole body including the head. This is a single patient use sling that tolerates water – dispose after discharge. Beneficial for infection control program. Same design as Guldmann ABC Basic High sling.

Suitable for

Lifting persons without head and body control • lifting to and from a lying position (bed, floor, etc) • lifting to and from a seated position.

Packed in boxes of 10 slings.

Material

Polyester – special woven

M/L 28366
L/XL 28367

2

Disposable Micro Plus

Function

Supports the upper part of the body (from upper pelvis till under the shoulders) around the chest and upper thigh. The belt has a re-tightening function if necessary to adjust the support around the body. This is a single patient use sling that tolerates water. Dispose after discharge. Same design as Guldmann ABC Active Micro Plus sling.

Suitable for

- lifting persons with head control, but with reduced muscular strength of upper part of the body and above hip joint and thighbone
- in situations where the sling should be put on and taken off easily
- lifting to and from seated position – in certain situations also lifting to and from a half-lying position together with Guldmann Head Support.

Packed in boxes of 10 slings.

Material

Polyester – special woven

XS 281535
S 281545
M 281555
L 281565
XL 281575
XXL 281585

Disposable slings



3
Disposable Leg sling

Used when legs need lifting and support. Possible activities are dressings, apply therapeutic sock/hose, access to hygiene, support of leg to apply ROM devices. Dispose after discharge. Beneficial for infection control program. Same design as Guldmann Leg sling.
Packed in boxes of 10 slings + 1 stiffener.

Material

Polyester – special woven

..... 286505

4
Disposable Repositioning sling, 455 kg

A sling for use as part of the linen make-up of the patient bed; leave positioned underneath patient, ready to use for care tasks. This is a single patient use sling that tolerates water – dispose after discharge. Beneficial for infection control program. Same design as Guldmann Repositioning sling.

Suitable for

Rolling to sidelying, repositioning patient up in bed, sitting forward head and trunk, supine lateral transfers to another surface such as beds + stretchers, and full lift in supine off of bed to zero bed scale, change linen, or other.
Packed in boxes of 10 slings.

Material

Polyester – special woven

..... 284655

Slings for 4 Point Hanger

Guldmann™

1



2



3



1

Vip Stretch

Sling for the 4 Point Hanger.
Used when lifting a person where it is required to change the seating position during the lifting process.

Material

Polyester

M/L 285151
L/XL 285161

2

Standard Net

Sling for the 4 Point Hanger.
Used when lifting a person where it is required to change the seating position during the lifting process, e.g. in connection with bathing.

Material

Net

M/L 285052
L/XL 285062

3

Vip Net

Sling for the 4 Point Hanger.
Used when lifting a person where it is required to change the seating position during the lifting process, e.g. in connection with bathing. This sling provides the user a more reclining position compared with the Standard Net sling.

Material

Net

M/L 285252
L/XL 285262




1
GH3 Lifting hanger, medium, 350 kg
 Suitable for persons using a lifting sling
 in sizes M – XXXL.
 550316



2
GH3 Lifting hanger, small, 350 kg
 Suitable for persons using a lifting sling in
 sizes XS and S.
 550057



3
GH3 Lifting hanger, x-small, 350 kg
 Suitable for persons using a lifting sling
 in sizes Kids and XS.
 550058

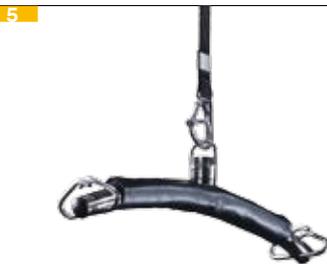

4
Lifting hanger, GH2, standard, 250 kg
 Suitable for persons using a lifting sling in sizes S-XXXL.
 for mobile lifters 21180
 for ceiling hoist w/fixed hanger 21170
 for ceiling hoist w/snap hook strap . . 21187

5
Lifting hanger, GH2, small, 250 kg
 Suitable for persons using a lifting sling in sizes XS and S.
 for mobile lifters 21280
 for ceiling hoist w/fixed hanger 21270
 for ceiling hoist w/snap hook strap . . 21287

6
Lifting hanger, GH2, x-small, 250 kg
 Suitable for persons using a lifting sling in sizes Kids and XS.
 for mobile lifters 21580
 for ceiling hoist w/fixed hanger 21570
 for ceiling hoist w/snap hook strap . . 21587

7
GH3 Cross hanger, 500 kg
 Suitable when larger spacing between back and leg straps of the lifting sling is required. To be used when lifting with GH3+ and GH3 Twin ceiling hoist.
 550444
 Connecting bar for GH3 Twin, 500 kg
 550544

8
H Hanger, 250 kg, GH2 + GH3
 Suitable when lifting persons where a large spacing between back and leg straps of the lifting sling is required. E.g. when lifting persons with a large abdomen/ stomach. The H hanger is available with side hangers in 3 different sizes.
 H Hanger – Centre hanger 21780
Note: In addition you need to order 2 hangers of either x-small (21587), small (21287) or standard (21187)

Lifting accessories

Guldmann™

9

Horizontal lifting support, 120 kg

Suitable when a horizontal lift is required.
With weight adjustment. To be used with
GH2, GH3 and mobile lifter.

..... 28450



10

Horizontal lifter, foldable, 200/250 kg

Suitable when a horizontal and stable lift
is required, with weight adjustment. To be
used with GH2 and GH3. Can be used with
either a complete lifting sheet or straps

..... 28445

Lifting sheet..... 284485

Set of lifting straps..... 284480



11


Horizontal lifter, foldable, 200/250 kg 2 strap version

Suitable when a horizontal and stable lift
is required, with weight adjustment. To be
used with GH2 HD Twin and GH3 Twin. Can
be used with either a complete lifting sheet
or straps

..... 28444

Lifting sheet..... 284485

Set of lifting straps..... 284480

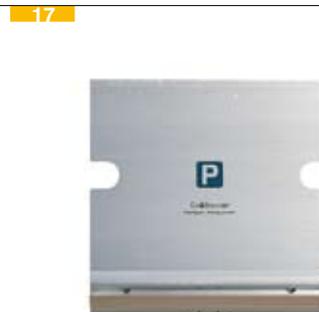

12
4 Point Hanger, 200 kg
 Suitable in situations where you wish/need to change the hanger's position during the lifting process. The hanger is available both in an electric and manual version.
 4 Point Hanger, electric for GH2 28500
 4 Point Hanger, manual for GH2 + GH3 and mobile lifter 28530
The hangers can only be used together with the slings for the 4-Point Hanger (p. 40).

13
Horizontal lifter BS, 150 kg
 A horizontal lifter especially developed and designed to overcome the many moving and handling challenges that mortuary staff faces every day. 28425

14
Bath chair – rustproof, 150 kg
 Suitable when a stable lift into a swimming pool is required. 14600

15
Bath stretcher – rustproof, 150 kg
 Suitable when a stable and horizontal lift into a swimming pool is required. 14610

16
Digital scale
 Battery operated digital scale.
 Max lifting capacity 320 kg.
 Approved for medical use.
 11346
 Max lifting capacity 250 kg 113461

17
Parking space for GH2 lifting hanger
 21140

18
Parking space for GH3 lifting hanger
 550531

19
Turning disc
 Allows a standing person to be rotated to correct position by eliminating the jerking movements that put stress on weak limbs and joints.
 50010

Lifting accessories



20
Leg sling
 Used when legs need lifting and support. Possible activities are dressings, apply therapeutic sock/hose, access to hygiene, support of leg to apply ROM devices.
Packed in boxes of 10 slings + 1 stiffener.
 28650

21
Turner
 Enables placing of sling under patient without having to roll or move them. Smaller surface area of material for this sling allows access to patient's skin for inspection, treatment and care 28700


22
Twin Turner
 Makes it possible to roll a patient on both sides without removing the sling. Enables placing of sling under patients without having to roll or move them. Smaller surface area of material for this sling allows access to patient's skin for inspection, treatment and care.
 Regular 28750
 Large 28751



23
Head support
 Supports the head without stimulating the neck reflexes. Suitable for persons with no head control. 28325

24
Sheepskin insets
 Suitable for persons who need more gentle support of the thighbone. 28315

25
Auxiliary straps
 Ease correct positioning of the user in the wheelchair. Use auxiliary straps to attach for handling. 28355

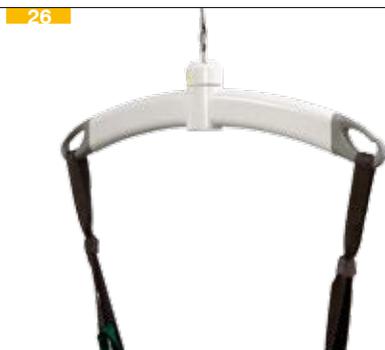

26
Extension straps
 Used in situations to achieve various sitting positions or where different lifting heights are required.
 20 cm 28470
 30 cm 28472

27
Transport trolley for GH2F
 12007



28
Repositioning sling, 455 kg
 A sling for use as part of the linen make-up of the patient bed; leave positioned underneath patient, ready to use for care tasks. Suitable for: rolling to sidelying, repositioning patient up in bed, sitting forward head and trunk, supine lateral transfers to another surface such as beds + stretchers, and full lift in supine off of bed to zero bed scale, stripe linen, or other. *Packed in boxes of 10 slings.* 28465



29

Trolley for horizontal lifter

..... 28443

taylordolman

tel: 0845 838 0600

email: info@taylordolman.com

